ahn my DP ette®

Bottle Top Dispenser

Unique piston mechanism - piston and cylinder can be cleaned without affecting the calibration

Dual graduation marks - volume can be read from both sides

Springless valve design - leads to smooth functioning and high chemical resistance

Eco-friendly - while discharging the tubes in the re-circulation mode, in fast priming, the reagent is redirected back into the bottle, eliminating any waste of precious reagents and ensuring safety and savings

Easy handling - discharge tube with 360° rotation, makes the bottle label visible all the time / telescopic filling tube compatible with different sizes of reagent bottles / ergonomic

Easy to calibrate - complies with ISO 8655-5 standard



Ordering information

Description	Sales Unit	Cat. No.
AHN myDpette Bottle Top Dispenser 0.25-2.5 mL	Carton/1 pc.	8-700-00-0
AHN myDpette Bottle Top Dispenser 0.5-5.0 mL	Carton/1 pc.	8-701-00-0
AHN myDpette Bottle Top Dispenser 1.0-10.0 mL	Carton/1 pc.	8-702-00-0
AHN myDpette Bottle Top Dispenser 2.5-25.0 mL	Carton/1 pc.	8-703-00-0
AHN myDpette Bottle Top Dispenser 5.0-50.0 mL	Carton/1 pc.	8-704-00-0


YEARS

autoclavable 120°C 20 min

Delivery Package
1 pc. Bottle top dispenser
4 pcs. Adapters (28, 38, 40 and 45mm screw cap)
1 pc. Telescopic filling tube
1 pc. Re-circulation tube
1 pc. Calibration tool
1 pc. User manual
1 pc. Certificate of conformity
1 pc. Calibration report
1 pc. Warranty card





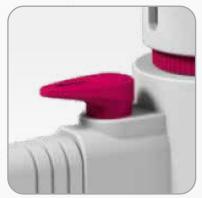

Recirculation valve











