

# AHN myLab® TR-03

Plate Rotator 30 rpm

**Usage** - ideal for processing of blood samples, separating and mixing of liquid and powdered samples

**Capacity** - 12 x 15 mL and 4 x 50 mL tubes or 22 x 15 mL tubes with the maximum operational weight capacity of 1.0 kgs

**Fixed speed** - 30 rpm with an accuracy of +/-2 rpm

**Optimal mixing** - circular motion allows gentle but effective mixing

**Robust and space-saving design** - strong metal construction with rubber feet for stable operation / compact design requires minimum of bench space

**Easy maintenance and adjustment** - plates can be easily removed and installed for cleaning and replacing



## Ordering information

### Description

AHN myLab® TR-03 Plate Rotator 30 rpm

### Cat. No.


7-220-00-0

### Delivery package

- 1 pc. Plate Rotator 30 rpm
- 1 pc. Power supply adaptor
- 1 pc. 12 x 15 mL + 4 x 50 mL Plate
- 1 pc. Instruction manual
- 1 pc. Warranty card

### Specifications

Run time	Continuous / infinite
Noise level	Noiseless operation
Max. loading capacity	1.0 kg
Protection class	IP 21
Dimensions (WxDxH)	346 x 185 x 145 mm
Weight	2.5 kg

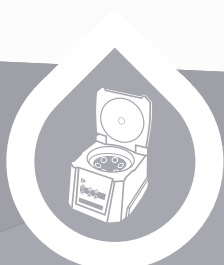

+49(0)3631/46594-04  
+49(0)3631/46594-10  
info@ahn-bio.de

Uthleber Weg 14  
D-99734 Nordhausen  
www.ahn-bio.de

ahn®  
biotechnologie GmbH