

ahn

2020/21
CATALOGUE



ahn

For over twenty years AHN Biotechnologie GmbH has been an OEM manufacturer of consumables for many top life science brands and distribution companies in Germany and worldwide. Located in Nordhausen, Germany, our state of the art production facility meets and exceeds the strict requirements of the most well-known brands we are supplying.

Now, it is time for AHN Biotechnologie GmbH to introduce its own brand and provide value directly to scientists worldwide. Besides our current products we are strengthening our brand by introducing a comprehensive range of consumables, liquid handling instruments and laboratory equipment. Our unique combination of flexibility, innovation and high level manufacturing automation will ensure laboratory products fulfilling the needs of worldwide customers and complying with the quality they desire.

Our know-how within pipette tips production has now been strengthened by the implementation of newly developed Swiss molds, leading to the launch of new tip-lines, matching the strict standards of modern lab industry. Our goal is to offer the markets best value for money on pipette tips and other AHN products.

To obtain the most optimal performance of our consumables, we have been continuously working on the improvement of our instruments. The new pipet4u® pro pipettes, with improved accuracy and greater user-friendliness along with the brand new AHN myDpette bottle top dispensers, convince by simplicity of use combined with the highest quality of workmanship. With our pipet4u® Revolution we offer an affordable range of electronic pipettes. Profit from the benefits of electronic dispensing, such as improved precision and multi dispensing function.

The AHN myLab® products are being brought to market with one thing in mind: YOU. A wide range of centrifuges, stirrers and shakers are available in a variety of sizes, multiple accessories and more features than most of the competitors.

We invite you to explore the new AHN and discover a better way to source your laboratory equipment and consumables.

www.ahn-bio.de

List of Contents



CONSUMABLES

Pipette Tips	4-5
Pipette Filter Tips	6-7
Low Retention Tips Filter Tips	8-9
Reload System	10
Rack Box System	11
Compatibility Chart	12-13
PCR Tubes	14
PCR Strips	15
PCR Plates	16
Sealing Membranes	17
Micro Tubes	18
Centrifuge Tubes	19
Reagent Reservoirs	20
Cryo Tubes	21
Cryoboxes	22-23



PURIFICATION

Spin Column System	24-25
Microfilterplate System	24-25



PIPETTING

Manual Pipettes	26-28
Pipette Stands	29
Electronic Pipettes	30-31
Pipette Controller	32
Bottle Top Dispenser	33



BENCHTOP EQUIPMENT

Microcentrifuges	34-35
Clinical Centrifuge	36
Hematocrit Centrifuge	37
Universal Centrifuge	38-39
Mini Vortex Mixer	40
Vortex Mixer	41
Stirrers	42-43
Tube Roller	44
Tube Rotators	45-46

ahn myTip

Pipette Tips | Macro Tips



4

ahn myTip Pipette Tips

Optimised fit and easy tip ejection - improved cone geometry for a secure fit

Quick and easy volume check - graduation marks on tip at standard volumes

Precise pipetting without loss of sample - conical tip opening prevents drop formation at the end of the tip


High quality material - tips made of autoclavable, non-wettable, virgin polypropylene

Easy identification of pipette tips - colour-coded racks allow a quick and easy determination of the pipette tip and the suitable pipette

Maintain sample integrity - sterile version available (beta irradiation)

Versions - available in bags and rack boxes

Compatibility list available on request

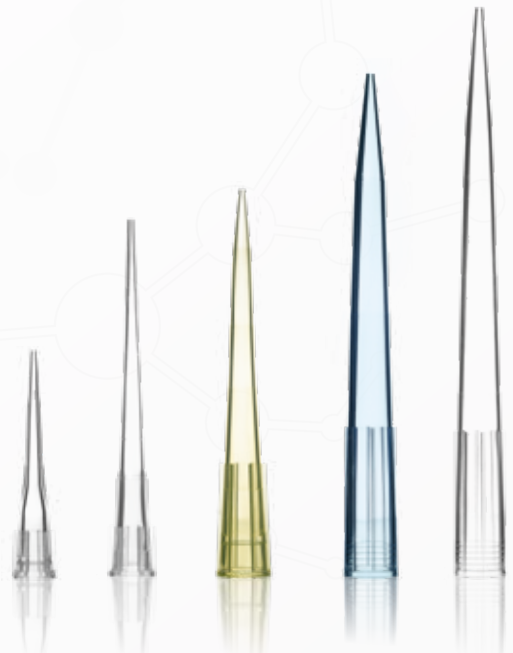

metal free



ahn myTip MaT Macro Tips

Optimised cone geometry - ensuring perfect fit on the defined pipette

Precision without loss of sample - conical tip opening to avoid drop formation

Elimination of sample contamination - tips made of autoclavable, non-wettable, virgin polypropylene material

Maintain sample integrity - sterile version available (beta irradiation)

Packaging - available in bags

Compatibility list available on request











metal free






ahn



Uthleber Weg 14 | 99734 Nordhausen | Germany

Ordering information

Volume	Length	Colour Tip	Colour Code	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® 10 µL	31.20 mm	clear			bag	case / 20 x 1000 tips	1-002-81-0
		clear	■		racked	case / 5 x 10 x 96 tips	1-002-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	1-003-C5-0
 AHN myTip® 10/20 µL, XL	45.80 mm	clear			bag	case / 20 x 1000 tips	1-062-81-0
		clear	■		racked	case / 5 x 10 x 96 tips	1-062-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	1-063-C5-0
 AHN myTip® 200 µL	53.40 mm	clear			bag	case / 30 x 1000 tips	1-121-86-0
		clear	■		racked	case / 5 x 10 x 96 tips	1-121-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	1-122-C5-0
 AHN myTip® 200 µL	53.40 mm	yellow			bag	case / 30 x 1000 tips	1-121-86-1
		yellow	■		racked	case / 5 x 10 x 96 tips	1-121-C5-1
		yellow	■	■	racked	case / 5 x 10 x 96 tips	1-122-C5-1
 AHN myTip® Gelloader 200 µL	68.50 mm	clear			bag	case / 5 x 1000 tips	1-130-70-0
		clear	■		racked	case / 4 x 8 x 96 tips	1-130-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	1-131-C4-0
 AHN myTip® 300 µL	59.35 mm	clear			bag	case / 10 x 1000 tips	1-127-80-0
		clear	■		racked	case / 5 x 10 x 96 tips	1-127-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	1-128-C5-0
 AHN myTip® 1000 µL, short	70.10 mm	clear blue			bag	case / 10 x 1000 tips	1-209-80-2
 AHN myTip® 1000 µL	84.30 mm	clear			bag	case / 10 x 1000 tips	1-202-80-0
		clear	■		racked	case / 4 x 8 x 96 tips	1-202-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	1-203-C4-0
 AHN myTip® 1000 µL	84.30 mm	blue			bag	case / 10 x 1000 tips	1-202-80-2
		blue	■		racked	case / 4 x 8 x 96 tips	1-202-C4-2
		blue	■	■	racked	case / 4 x 8 x 96 tips	1-203-C4-2
 AHN myTip® 1250 µL	97.60 mm	clear			bag	case / 10 x 1000 tips	1-204-80-0
		clear	■		racked	case / 4 x 8 x 96 tips	1-204-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	1-205-C4-0

Volume	Length	Colour	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® MaT 5 mL	124.00 mm	clear		bag	case / 10 x 250 tips	1-218-61-0
		clear	■	bag	case / 10 x 250 tips	1-222-61-0
 AHN myTip® MaT 5 mL	126.00 mm	clear		bag	case / 8 x 250 tips	1-210-84-0
 AHN myTip® MaT 5.5 mL	147.50 mm	clear		bag	case / 8 x 250 tips	1-212-84-0
		clear	■	bag	case / 8 x 250 tips	1-216-84-0
 AHN myTip® MaT 10 mL	150.00 mm	clear		bag	case / 10 x 100 tips	1-213-50-0
 AHN myTip® MaT 10 mL	181.00 mm	clear		bag	case / 5 x 200 tips	1-220-50-0
		clear	■	bag	case / 5 x 200 tips	1-224-50-0



ahn myTip FT

Pipette Filter Tips



metal free



6

No sample contamination - filters are made of virgin high quality hydrophobic polyethylene without additives

No pipette contamination - filter provides reliable protection against contamination through micropore filter technology (MPF technology, filter pore size 10-40 μm)

Universal fit and tip ejection - optimised cone geometry for a secure fit

Quick and easy volume check - graduation marks on tip at standard volumes

Precise pipetting without loss of sample - conical tip opening prevents drop formation at the end of the tip

Elimination of sample contamination - tips made of autoclavable, non-wettable, virgin polypropylene







Easy identification of pipette tips - colour-coded racks allow a quick and easy determination of the pipette tip and the suitable pipette

Maintain sample integrity - sterile version available (beta irradiation)









Versions - available in bags and rack boxes

Compatibility list available on request



Filter	Average pore size	Material	Characteristic	Used in following filter tip
	10-40 μm	PE	hydrophobic	10 μL 10 μL XL
	10-40 μm	PE	hydrophobic	20 μL
	10-40 μm	PE	hydrophobic	100 μL
	10-40 μm	PE	hydrophobic	200 μL 300 μL
	10-40 μm	PE	hydrophobic	1000 μL
	20-40 μm	PE	hydrophobic	1250 μL

Ordering information

Volume	Length	Colour Tip	Colour Code	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® FT 10 µL	31.20 mm	clear			bag	case / 20 x 1000 tips	2-002-81-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-003-C5-0
 AHN myTip® FT 10/20 µL XL	45.80 mm	clear			bag	case / 20 x 1000 tips	2-062-81-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-063-C5-0
 AHN myTip® FT 20 µL	53.40 mm	clear			bag	case / 10 x 1000 tips	2-117-80-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-118-C5-0
 AHN myTip® FT 100 µL	53.40 mm	clear			bag	case / 10 x 1000 tips	2-119-80-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-120-C5-0
 AHN myTip® FT 200 µL	53.40 mm	clear			bag	case / 10 x 1000 tips	2-127-80-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-128-C5-0
 AHN myTip® FT 300 µL	59.35 mm	clear			bag	case / 10 x 1000 tips	2-134-80-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	2-135-C5-0
 AHN myTip® FT 1000 µL	84.30 mm	clear			bag	case / 10 x 1000 tips	2-202-80-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	2-203-C4-0
 AHN myTip® FT 1250 µL	97.60 mm	clear			bag	case / 10 x 1000 tips	2-204-80-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	2-205-C4-0



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.dhn-bio.de



ahn myTip LT | LFT



Low Retention Tips | Filter Tips

8

Low retention - the innovative low retention surface ensures maximum sample recovery and high-precision pipetting

Universal fit and tip ejection - optimized cone geometry for a secure fit

Quick and easy volume check - graduation marks on tip at standard volumes

Precise pipetting without loss of sample - conical tip opening prevents drop formation at the end of the tip

Elimination of sample contamination - tips made of autoclavable, non-wettable, virgin polypropylene

Maintain sample integrity - sterile version available (beta irradiation)

No sample contamination - filters are made of virgin high quality hydrophobic polyethylene without additives

No pipette contamination - filter provides reliable protection against contamination through micropore filter technology (MPF technology, filter pore size 10-40 μm)

Easy identification of pipette tips - colour-coded racks allow a quick and easy determination of the pipette tip and the suitable pipette

Versions - available in bags and rack boxes

Compatibility list available on request








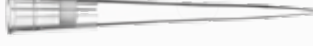



Standard Tip



AHN myTip® LT | LFT
Low Retention Tips | Filter Tips

Ordering information

Volume	Length	Colour Tip	Colour Code	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® LT 10 µL	31.20 mm	clear			bag	case / 20 x 1000 tips	4-002-81-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	4-003-C5-0
 AHN myTip® LT 10/20µL XL	45.80 mm	clear			bag	case / 20 x 1000 tips	4-062-81-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	4-063-C5-0
 AHN myTip® LT 200 µL	53.40 mm	clear			bag	case / 30 x 1000 tips	4-121-86-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	4-122-C5-0
 AHN myTip® LT 300 µL	59.35 mm	clear			bag	case / 10 x 1000 tips	4-127-80-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	4-128-C5-0
 AHN myTip® LT 1000 µL	84.30 mm	clear			bag	case / 10 x 1000 tips	4-202-80-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	4-203-C4-0
 AHN myTip® LT 1250 µL	97.60 mm	clear			bag	case / 10 x 1000 tips	4-204-80-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	4-205-C4-0

Volume	Length	Colour Tip	Colour Code	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® LFT 10 µL	31.20 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-003-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	5-063-C5-0
 AHN myTip® LFT 10/20 µL XL	45.80 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-118-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	5-120-C5-0
 AHN myTip® LFT 20 µL	53.40 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-128-C5-0
		clear	■	■	racked	case / 5 x 10 x 96 tips	5-135-C5-0
 AHN myTip® LFT 100 µL	53.40 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-203-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	5-205-C4-0
 AHN myTip® LFT 200 µL	59.35 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-203-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	5-205-C4-0
 AHN myTip® LFT 300 µL	84.30 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-203-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	5-205-C4-0
 AHN myTip® LFT 1000 µL	97.60 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-203-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	5-205-C4-0
 AHN myTip® LFT 1250 µL	97.60 mm	clear	■	■	racked	case / 5 x 10 x 96 tips	5-203-C4-0
		clear	■	■	racked	case / 4 x 8 x 96 tips	5-205-C4-0



ahn myTip[®] RL / RB

Reload System (RL)



10

Eco-friendly and cost-effective tip management solution – each tower accommodates multiple insert plates of refill tips, minimising plastic waste **by approximately 70% per pack** and reusing the current AHN myTip[®] RB rack boxes by simply exchanging empty plates with full ones from the AHN myTip[®] RL reload towers

Secure fit of pipette tips – pipette tips sit stable and secure in the insert plate, enabling easy handling with multichannel pipettes

Secure and comfortable handling – AHN myTip[®] RL insert plate securely snaps into the rackbox, but can be easily removed again. To facilitate the handling for users with small hands, the tabs of the AHN myTip[®] RL insert plate are located on the longer side

Easy identification of pipette tips – colour-coded racks allow quick and easy identification of the pipette tip, and therefore the suitable pipette

Space-saving – small footprint of the AHN myTip[®] RL reload towers means there is more space on the benchtop and in the storage room



metal free



1

Open the refill by sliding it towards the top

2



Prepare an empty rack. Remove its original grid

3

Pick up the first grid of refill tips and insert into the rack. The racks can be reused

4

Once the tips are properly installed in the rack, you may close the hinged lid

Volume	Length	Colour Tip	Colour Code	Sterilised	Pack Type	Sales Unit	Cat. No.
 AHN myTip® RL 10 µL	31.20 mm	Clear	■		reload	case / 10 x 10 x 96 tips	1-002-C1-R
 AHN myTip® RL 10/20 µL XL	45.80 mm	Clear	■		reload	case / 10 x 10 x 96 tips	1-062-C1-R
 AHN myTip® RL 200 µL	53.40 mm	Clear	■		reload	case / 10 x 10 x 96 tips	1-121-C1-R
 AHN myTip® RL 300 µL	59.35 mm	Clear	■		reload	case / 10 x 10 x 96 tips	1-127-C2-R
 AHN myTip® RL 1000 µL	84.30 mm	Clear	■		reload	case / 10 x 10 x 96 tips	1-202-C3-R

Rack Box System without Tips (RB)

Effortless and secure handling - easy one-handed opening and closing of rack box / secure rack box closing / rack boxes can be stacked for easy storage

Universal fit - suitable for multichannel pipettes (standard 9 mm SBS footprint)

Easy identification of pipette tip - colour coded insert plates allow for quick and easy identification of the pipette tip volume

Rack boxes available in small and large size / racks in 8 x 12 format



Ordering information

Rack Box System without Tips (RB)

Volume		Colour	Sales Unit	Cat. No.
AHN myTip® RB for tips 10 µL/10 µL XL		clear black (RB) / light grey (I)*	case / 5 x 10 rack boxes	4-500-C5-E
AHN myTip® RB for tips 200 µL/300 µL		clear black (RB) / yellow (I)*	case / 5 x 10 rack boxes	4-501-C5-E
AHN myTip® RB for tips 200 µL Gelloader		clear black (RB) / yellow (I)*	case / 4 x 8 rack boxes	4-502-C4-E
AHN myTip® RB for tips 1000 µL		clear black (RB) / blue (I)*	case / 4 x 8 rack boxes	4-503-C4-E
AHN myTip® RB for tips 1250 µL		clear black (RB) / blue (I)*	case / 4 x 8 rack boxes	4-504-C4-E

*RB - Rack Box | I - Insert Plate



Compatibility Chart

Description	Cat. No.	Gilson		Eppendorf		Brand		Thermo		Sartorius		HTL	
		Pipetman G single and multi channel pipettes	Research Plus single and multi channel pipettes	Transferring single and multi channel pipettes	Finnpipette F1 / Novus manual single and multi channel and electronic single and multi channel pipettes	Proline Plus / Picuss manual single channel and electronic single and multi channel pipettes	Labmate Pro (same as Labnet BioPette) manual single and multi channel pipettes						
AHN myTip® Pipette Tips													
AHN myTip® 10 µL, clear	1-002-81-0	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL	10 µL
AHN myTip® 10 µL, clear	1-002-C5-0	1000-10000 µL	1000-10000 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 10 µL, clear, sterile	1-003-C5-0	1000-5000 µL	1000-5000 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 10/20 µL, clear	1-042-81-0	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL	20 µL
AHN myTip® 10/20 µL, clear	1-042-C5-0	100-1000 µL	100-1000 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 10/20 µL, clear, sterile	1-043-C5-0	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (grey)	2-20 µL (grey)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)	2-20 µL (yellow)
AHN myTip® 200 µL, clear	1-121-86-0	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL	0.2-2 µL
AHN myTip® 200 µL, clear	1-121-C5-0	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL	300 µL
AHN myTip® 200 µL, clear, sterile	1-122-C5-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 200 µL, yellow	1-121-86-1	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 200 µL, yellow	1-121-C5-1	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 200 µL, yellow, sterile	1-122-C5-1	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 200 µL, Gelloader, clear	1-130-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 200 µL, Gelloader, clear, sterile	1-131-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 300 µL, clear	1-127-80-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 300 µL, clear	1-127-C5-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 300 µL, clear, sterile	1-128-C5-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, short, clear blue	1-209-80-2	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, clear	1-202-80-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, clear	1-202-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, clear, sterile	1-203-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, blue	1-202-80-2	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, blue	1-202-C4-2	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1000 µL, blue, sterile	1-203-C4-2	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1250 µL, clear	1-204-80-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1250 µL, clear	1-204-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® 1250 µL, clear, sterile	1-205-C4-0	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
AHN myTip® Mat Macro Pipette Tips													
AHN myTip® Mat 5 mL, clear	1-218-61-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 5 mL, clear, sterile	1-222-61-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 5 mL, clear	1-210-84-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 5.5 mL, clear	1-212-84-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 5.5 mL, clear, sterile	1-216-84-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 10 mL, clear	1-213-50-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 10 mL, clear	1-220-50-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL
AHN myTip® Mat 10 mL, clear, sterile	1-224-50-0	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL	1000-10000 µL

■ compatible x not compatible

The compatibility chart above is for your reference. We are not responsible for technical modifications or changes made to the mentioned pipettes. Specifications like inaccuracy and imprecision values have not been tested. For additional compatibility charts, please contact info@cappahn.com.



Compatibility Chart



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



Description	Cat. No.	Gilson		Eppendorf		Brand		Thermo						Sartorius				HTL																						
		Pipetman / Pipetman G	Research Plus	Transferpipette, Transferpipette S	single and multi channel pipettes	single and multi channel pipettes	single and multi channel pipettes	Finnpipette F1 / Novus	manual single and multi channel pipettes	electronic multi	electronic	single	electronic	electronic multi	manual single	manual single	manual multi	Labmate Pro (same as Labnet BioPipette)	manual single and multi channel pipettes																					
AHN myTip® Filter Tips		0.2-2 µl	2-20 µl	20 µl	300 µl	10 µl	100 µl	1000-10000 µl	1000-10000 µl	100 µl	10 µl	0.5-10 µl pink	0.5-10 µl yellow	2-20 µl yellow	2-20 µl turquoise	10-100 µl	20-200 µl	300 µl	0.5-10 µl	20-200 µl	300 µl	100-1000 µl	100 µl	10 µl	0.5-10 µl	1000-10000 µl	500-5000 µl	1000-10000 µl	50-300 µl											
AHN myTip® FT 10 µL, clear	2-002-81-0		X									X																												
AHN myTip® FT 10 µL, clear, sterile	2-003-C5-0		X									X																												
AHN myTip® FT 10/20 µL, clear	2-042-81-0		X									X																												
AHN myTip® FT 10/20 µL, clear, sterile	2-043-C5-0		X									X																												
AHN myTip® FT 20 µL, clear	2-117-80-0		X									X																												
AHN myTip® FT 20 µL, clear, sterile	2-118-C5-0		X									X																												
AHN myTip® FT 100 µL, clear	2-119-80-0		X									X																												
AHN myTip® FT 100 µL, clear, sterile	2-120-C5-0		X									X																												
AHN myTip® FT 200 µL, clear	2-127-80-0		X									X																												
AHN myTip® FT 200 µL, clear, sterile	2-128-C5-0		X									X																												
AHN myTip® FT 300 µL, clear	2-134-80-0		X									X																												
AHN myTip® FT 300 µL, clear, sterile	2-135-C5-0		X									X																												
AHN myTip® FT 1000 µL, clear	2-202-80-0		X									X																												
AHN myTip® FT 1000 µL, clear, sterile	2-203-C4-0		X									X																												
AHN myTip® FT 1250 µL, clear	2-204-80-0		X									X																												
AHN myTip® FT 1250 µL, clear, sterile	2-205-C4-0		X									X																												
AHN myTip® LT Low Retention Pipette Tips																																								
AHN myTip® LT 10 µL, clear	4-002-81-0		X									X																												
AHN myTip® LT 10 µL, clear, sterile	4-003-C5-0		X									X																												
AHN myTip® LT 10/20 µL, clear	4-062-81-0		X									X																												
AHN myTip® LT 10/20 µL, clear, sterile	4-063-C5-0		X									X																												
AHN myTip® LT 200 µL, clear	4-121-86-0		X									X																												
AHN myTip® LT 200 µL, clear, sterile	4-122-C5-0		X									X																												
AHN myTip® LT 300 µL, clear	4-127-80-0		X									X																												
AHN myTip® LT 300 µL, clear, sterile	4-128-C5-0		X									X																												
AHN myTip® LT 1000 µL, clear	4-202-80-0		X									X																												
AHN myTip® LT 1000 µL, clear, sterile	4-203-C4-0		X									X																												
AHN myTip® LT 1250 µL, clear	4-204-80-0		X									X																												
AHN myTip® LT 1250 µL, clear, sterile	4-205-C4-0		X									X																												
AHN myTip® LFT Low Retention Filter Tips																																								
AHN myTip® LFT 10 µL, clear, sterile	2-003-C5-0		X									X																												
AHN myTip® LFT 10/20 µL, clear, sterile	2-043-C5-0		X									X																												
AHN myTip® LFT 200 µL, clear, sterile	2-118-C5-0		X									X																												
AHN myTip® LFT 100 µL, clear, sterile	2-120-C5-0		X									X																												
AHN myTip® LFT 200 µL, clear, sterile	2-128-C5-0		X									X																												
AHN myTip® LFT 300 µL, clear, sterile	2-135-C5-0		X									X																												
AHN myTip® LFT 1000 µL, clear, sterile	2-203-C4-0		X									X																												
AHN myTip® LFT 1250 µL, clear, sterile	2-205-C4-0		X									X																												

■ compatible x not compatible

The compatibility chart above is for your reference. We are not responsible for technical modifications or changes made to the mentioned pipettes. Specifications like inaccuracy and imprecision values have not been tested. For additional compatibility charts, please contact info@cappahn.com.

ahn myTube PCR



PCR Tubes 0.1 mL | 0.2 mL

14

Elimination of contamination and loss of sample - high quality raw material meets purity standards, is extremely resistant to extreme mechanical strain and has excellent chemical resistance

Improved sample recovery - the V-shape bottom guarantees easier sample recovery

Reduced cycle time and optimal heat transfer - guaranteed by the extra thin, uniform wall thickness of the tube

Low-binding - the 0.1 mL PCR tube has a low-binding surface that reduces the adhesion of proteins to the very minimum

Optical measurement - high transparency of PCR tubes allows the use for applications where the optical qPCR signal comes from the bottom of the tubes

Centrifuge specifications: up to 18.000 x g - fixed angle rotor up to 500 x g - swing out rotor

Effortless and secure working - ergonomic flat cap design allows opening and closing with one hand

Sterile handling - each lot is manufactured under strict sterile conditions

Versions - available in two sizes - 0.1 mL and 0.2 mL



V-shaped bottom



Flat cap design

Ordering information

Volume	Height	Colour	Sterilised	Pack Type	Sales Unit	Cat. No.
0.1 mL 	16.6 mm	clear		bag	20x1000	3-000-80-0
0.2 mL 	21.5 mm	clear		bag	10x1000	3-010-80-0
		clear	■	bag	10x1000	3-011-80-0

metal free

autoclavable
121°C
15 min

ahn myTube PCR



PCR and qPCR Strips 8 x 0.1 mL | 8 x 0.2 mL

Optical measurement - high transparency of PCR strip (3-036-01-0) allows the use for applications where the optical qPCR signal comes from the bottom of the tubes / ultra clear cap shows consistent light transmitting rate

Compatibility - the 0.1 mL strips are especially recommended for fast thermocyclers requiring low profile tubes and the 0.2 mL strips for the standard thermocyclers requiring regular profile tubes / compatibility chart available on request

Reduced cycle time and optimal heat transfer - guaranteed by extra thin, uniform wall thickness of tube

Low-binding - 3-036-02-0 has a low binding surface that reduces the adhesion of proteins to the very minimum

Sterile handling - each lot is manufactured under strict sterile conditions

Reproducible result - evaporation rate below 5%

Versions - available as 0.1 and 0.2 mL PCR strip tubes, with separate caps or with attached caps



metal free

autoclavable
121°C
15 min

Ordering information

Volume		Height	Colour	Pack Type	Sales Unit	Cat. No.
8 x 0.1 mL (qPCR) with separate optical flat caps		16.0 mm	clear	bag	5x10x12	3-036-01-0
8 x 0.1 mL (PCR and qPCR) with attached caps		16.6 mm	clear	bag	8x120	3-036-02-0
8 x 0.1 mL (qPCR) with separate optical flat caps		16.0 mm	white	bag	5x10x12	3-036-03-9
8 x 0.2 mL (PCR) with separate frosted flat caps		21.0 mm	clear	bag	5x10x12	3-037-01-0
8 x 0.2 mL (PCR) with attached frosted flat caps		21.0 mm	clear	bag	5x10x12	3-037-02-0
8 x 0.2 mL (qPCR) with separate optical flat caps		21.0 mm	clear	bag	5x10x12	3-037-03-0
8 x 0.2 mL (qPCR) with separate optical flat caps		21.0 mm	white	bag	5x10x12	3-037-04-9



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



ahn myPlate PCR



PCR and qPCR Plates 96 x 0.1 mL | 96 x 0.2 mL

16

Optimal format - fits in various thermocyclers

Optical measurement - high transparency of PCR plate (3-046-00-0) allows the use for applications where the optical qPCR signal comes from the bottom of the plate

Improved sample recovery - the V-shape bottom guarantees easier sample recovery

Low-binding - 3-046-00-0 has a low binding surface, that reduces the adhesion of proteins to the very minimum


Sterile handling - each lot is manufactured under strict sterile conditions

Sealable - using disposable adhesive for secure closure of PCR plate (see page 17)

Versions - available as low profile (0.1 mL) and as regular profile (0.2 mL) PCR plates and as non-skirted, semi-skirted and sub semi-skirted format

Reproducible result - evaporation rate below 5%

Compatibility chart available on request

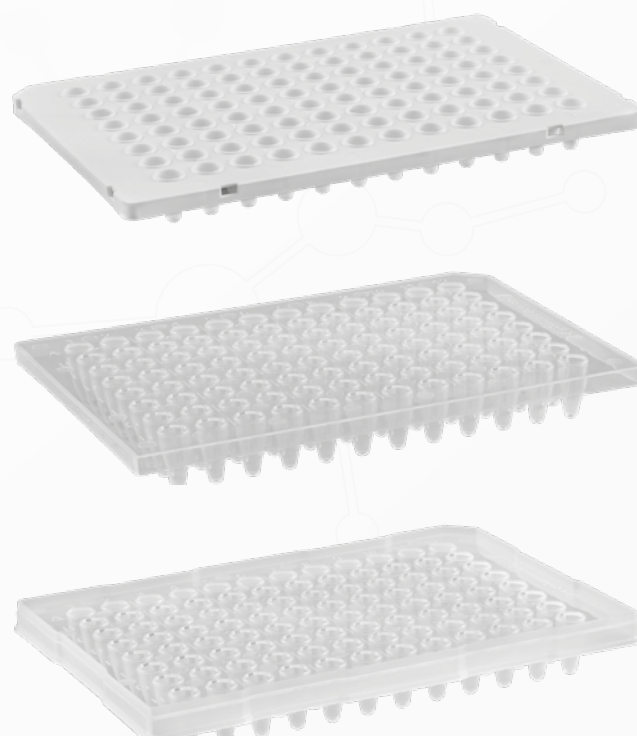

metal free

autoclavable
121°C
15 min

Ordering information

Volume		Height	Colour	Pack Type	Sales Unit	Cat. No.
96 x 0.1 mL (PCR and qPCR) non-skirted		16.6 mm	clear	bag	8x25	3-046-00-0
96 x 0.1 mL (qPCR) semi-skirted		16.0 mm	frosted	bag	5x5x10	3-046-01-0
96 x 0.1 mL (qPCR) sub semi-skirted		16.0 mm	frosted	bag	5x5x10	3-046-02-0
96 x 0.1 mL (qPCR) semi-skirted		16.0 mm	white	bag	5x5x10	3-046-03-9
96 x 0.2 mL (PCR) non-skirted		21.0 mm	clear	bag	5x5x10	3-047-01-0
96 x 0.2 mL (qPCR) semi-skirted		21.0 mm	clear	bag	5x5x10	3-047-02-0
96 x 0.2 mL (qPCR) sub semi-skirted		21.0 mm	clear	bag	5x5x10	3-047-03-0

ahn myPlate PCR




PCR and qPCR Sealing Membranes

Usage - designed for PCR and qPCR applications and for storage

Reproducible result - evaporation rate below 5 %


Ultra clear membranes - the optical sealing membranes 3-048-01-0 and 3-048-02-0 show consistent and high light transmitting rate for stable data

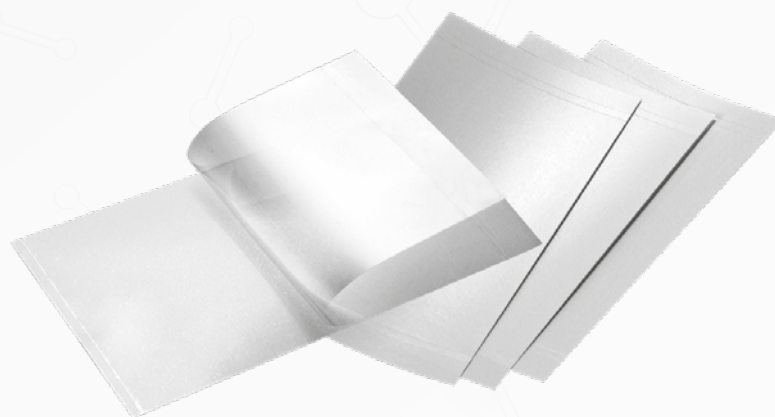

17

Ordering information

Volume	Features	Pack Type	Sales Unit	Cat. No.
96 Well Plate Sealing Membrane 	<ul style="list-style-type: none"> operating temperature -40 ~120°C PCR use strong adhesive 	bag	10x10	3-048-00-0
96 Well Plate Optical Sealing Membrane 	<ul style="list-style-type: none"> operating temperature -40 ~120°C Real-Time PCR(qPCR) use strong pressure sensitive adhesive 	bag	10x10	3-048-01-0
96 Well Plate Optical Sealing Membrane 	<ul style="list-style-type: none"> operating temperature -40 ~120°C Real-Time PCR(qPCR) use strong adhesive 	bag	10x10	3-048-02-0
Aluminium Membrane 	<ul style="list-style-type: none"> operating temperature -80 ~120°C PCR or Storage use strong adhesive piercable 	bag	10x10	3-048-03-0

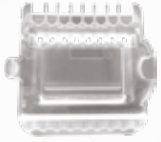Accessory	Features	Pack Type	Sales Unit	Cat. No.
Decapping/ Sealing Tool for PCR Strips PCR Membranes 	<ul style="list-style-type: none"> PCR strip decapping: 3-037-01-0, 3-037-02-0, 3-037-03-0, 3-037-04-9 Membrane sealing: 3-048-00-0, 3-048-01-0, 3-048-02-0, 3-048-03-0 	bag	15x1	3-000-OP-0



ahn myTube MCT



Micro Tubes 0.5 mL | 1.5 mL | 2.0 mL

18

Effortless and secure handling - ergonomic cap design with hinge orientation nibs (1.5 mL / 2.0 mL) and safe click lock, allowing an easy and secure opening and closing of the tube / bottom and cap have no sharp points, gloves will not be damaged / cap of the tube is puncturable

Quick and simple volume check - tubes are equipped with precise graduation marks / 1.5 mL tube with 0.2 mL extra space and 0.5 mL tube with 0.1 mL extra space for freezing / boiling purposes

Easy sample identification - frosted marking area for easy labeling

Optical measurement and unsurpassed clarity

Elimination of contamination and loss of sample - high quality raw material meets standards of purity and is highly resistant to mechanical strain and has an excellent chemical resistance

Centrifuge specifications: up to 20.000 x g - 0.5 mL tube up to 25.000 x g - 1.5 mL / 2.0 mL tube

Maintain sample integrity - sterile version available

Sterile handling - micro tubes are autoclavable (opened 121°C / 15 min)






metal free

CE IVD

autoclavable
121°C
15 min

Ordering information

Volume	Height	Colour	Sterilised	Pack Type	Sales Unit	Cat. No.
0.5 mL 	31.6 mm	clear	■	bag	5x1000	3-116-C5-0
		clear	■	bag	16x500	3-205-80-0
1.5 mL 	40.9 mm	clear	■	bag	5x1000	3-207-C5-0
		black*	■	bag	16x500	3-205-80-7
2.0 mL 	41.7 mm	clear	■	bag	16x500	3-204-80-0
		clear	■	bag	5x1000	3-208-C5-0

*protection of light sensitive samples

ahn myTube CT

Centrifuge Tubes 15 mL | 50 mL



Elimination of contamination of sample - high quality raw material meets purity standards

Centrifuge Tube with screw flat cap

Easy sample identification - large writing area on tube

Quick and simple volume check - precise graduation marks

No loss of sample - tube resists high centrifugal force and temperature variations (-80°C to +121°C)

Centrifuge specifications:

up to 17.000 x g - 15 mL centrifuge tube

up to 20.000 x g - 50 mL centrifuge tube

up to 7.000 x g - 50 mL centrifuge tube (self-standing)

Comfortable handling - the 3/4th turn closure enables one-handed operation by opening and tightening

Maintain sample integrity - sterile version available (by gamma irradiation)

Sterile handling - centrifuge tubes are autoclavable (opened 121°C / 15 min)

Versions - 50 mL tubes are available as non-skirted or self-standing (skirted)



Ordering information

Volume	Height	Colour	Sterilised	Pack Type	Sales Unit	Cat. No.
15 mL	120.0 mm	clear/grey		bag	20x25	3-254-25-8
		clear/grey		rack	20x25	3-255-25-8
		clear/grey	■	bag	20x25	3-256-25-8
		clear/grey	■	rack	20x25	3-257-25-8
50 mL	115.0 mm	clear/grey		bag	20x25	3-270-25-8
		clear/grey		rack	20x25	3-271-25-8
		clear/grey	■	bag	20x25	3-272-25-8
		clear/grey	■	rack	20x25	3-273-25-8
50 mL self-standing	115.0 mm	clear/grey		bag	20x25	3-280-25-8
		clear/grey	■	bag	20x25	3-281-25-8



Non-skirted



Self-standing (skirted)



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



ahn Reagent Reservoirs

10 mL | 25 mL | 50 mL

20

Construction - made out of a unique blend of Polystyrene

Extra large, sturdy base with broad edges - enable it to stand firmly on your workbench

Graduation marks - for easy, accurate dispensing / helps identify any remaining liquid

Angled interior supports the collection of the liquid and decreases reagent waste

Spillproof - no need to worry about spillage while you are processing liquids

Easy pouring - to store any remaining liquid you can empty the reservoir using the corners

Compact storage - easy to stack the reservoirs and separate them using the reservoir lip



Ideal for using with all AHN Single and Multichannel Pipettes

- to suit the various pipettes our reservoirs come in three different sizes:

- ① 10 mL, V-shaped bottom (ideal with 8- and 16-channel pipettes)
- ② 25 mL, flat bottom (ideal with multichannel pipettes)
- ③ 50 mL, V-shaped bottom (ideal with 12-channel pipettes)



Ordering information

Description	Pack Type	Sales Unit	Cat. No.
AHN Reagent Reservoirs 10 mL 	bag	5x100	3-701-11-9
AHN Reagent Reservoirs 25 mL 	bag	5x100	3-702-11-9
AHN Reagent Reservoirs 50 mL 	bag	5x100	3-703-11-9

ahn myTube CrT

Cryo Tubes 0.5 mL | 1.0 mL | 2.0 mL | 4.0 mL | 5.0 mL

Super Straight design - external thread closure saves up to 23% storage capacity

Quick Seal design - optimized for single-handed operation. Open or close with only half a turn

Leak-proof cap - the lip seal cap design ensures that the tubes are leak-proof

Quick and easy volume check - graduation marks on cryo tube with 0.5 mL intervals (not available for 0.5 mL cryo tube)

Easy sample identification - big writing area for easy sample identification (not available for 0.5 mL cryo tube)

Better stability and operation - the star-foot bottom improves the stability and operation of our cryo tubes when they are stored in selected boxes

autoclavable
121°C
15 min



Ordering information

Volume	Height	Colour	Sterilised	Pack Type	Sales Unit	Cat. No.
0.5 mL	28.6 mm	clear	■	bag	10x100	3-242-51-0
1.0 mL	36.1 mm	clear	■	bag	10x100	3-243-51-0
2.0 mL	47.6 mm	clear	■	bag	10x100	3-244-51-0
2.0 mL	47.6 mm	clear	■	bag	2x500 tubes / 1x1000 caps	3-244-50-0
4.0 mL	73.3 mm	clear	■	bag	5x100	3-245-51-0
5.0 mL	88.5 mm	clear	■	bag	5x100	3-246-51-0



ahn Cryoboxes

Plastic Cryoboxes 9x9 | 10x10

22

Store samples effortlessly – organise your cryo tubes in freezers or liquid nitrogen

Construction – plastic cryoboxes made from premium quality polypropylene

Broad thermal range – easy to autoclave

Two sizes available – 9x9 and 10x10 grid

Standard footprint – 130x130 mm, stacks evenly. Fits in regular freezers


Streamlined design – complements the AHN myTube CrT Cryo Tubes and retains compatibility with other brands

9x9 Cryoboxes – accommodate up to 81 samples. Suited to budget-constrained laboratory settings

10x10 Cryoboxes – provide 23% more storage without increasing the box footprint. Suited to space-limited settings or laboratories requiring cutting-edge features for optimizing workflow


For more colour options please contact info@cappahn.com

Additional colour options 

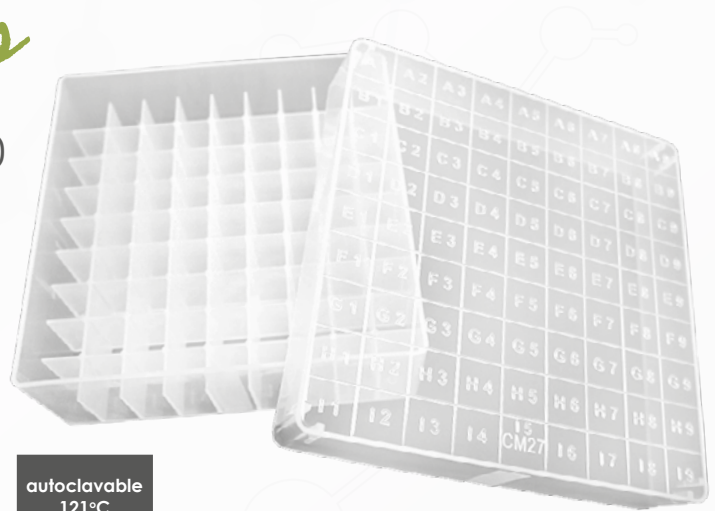

autoclavable
121°C
15 min

Features	9x9 Cryoboxes	10x10 Cryoboxes
Grid	9x9 square dividers	10x10 square dividers (with tube lock system)
Capacity	81 cryogenic tubes	100 cryogenic tubes
Material	Premium quality polypropylene	Premium quality polypropylene
Footprint	130x130 mm	130x130 mm
Box height	50, 80 or 90 mm	32, 40, 50, 80 or 90 mm
Colour	Clear	Clear
Thermal range	-90° to +121°C (autoclavable)	-150° to +121°C (autoclavable)
Sample sorting	Alpha-numeric system on top and bottom	Alpha-numeric system on top and bottom
Extras		One handed opening and closing of sample tubes. Slanted edges for quick, accurate closing

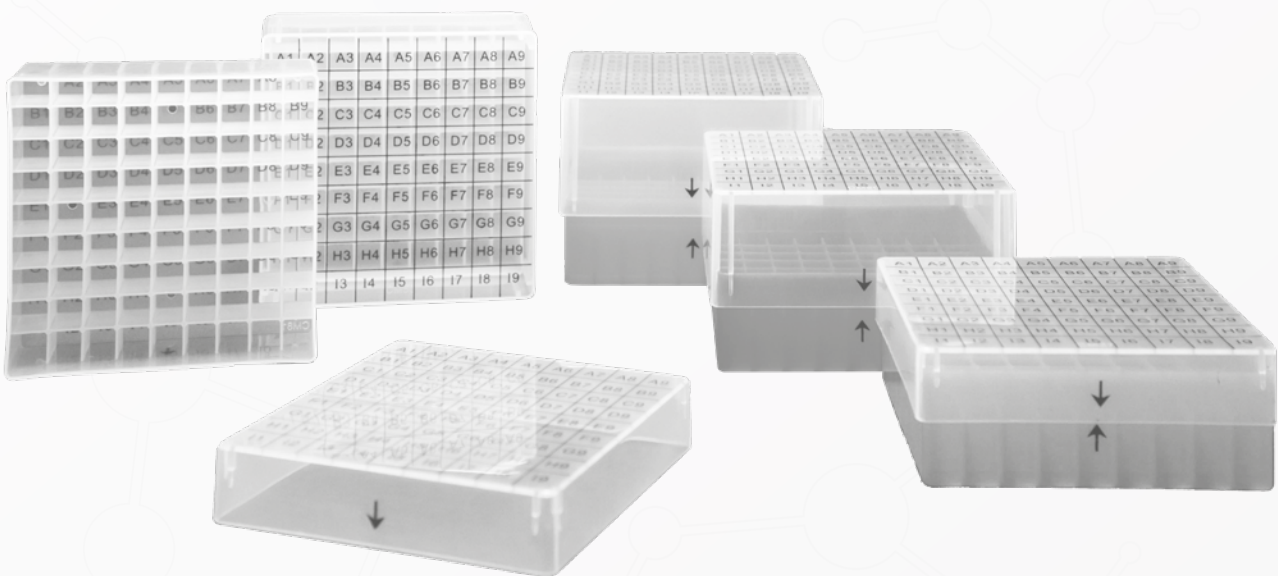

Description	Dimensions	Material	Colour	Sales Unit	Cat. No.
Cryobox, 9 x 9 grid for 1.0 - 2.0 mL cryo tubes	130x130x50 mm	Plastic	clear	case / 36 pcs.	3-632-45-0
Cryobox, 9 x 9 grid for 4.0 mL cryo tubes	130x130x80 mm	Plastic	clear	case / 36 pcs.	3-633-36-0
Cryobox, 9 x 9 grid for 5.0 mL cryo tubes	130x130x90 mm	Plastic	clear	case / 36 pcs.	3-634-36-0
Cryobox, 10 x 10 grid for 0.5 mL cryo tubes	130x130x32 mm	Plastic	clear	case / 36 pcs.	3-635-45-0
Cryobox, 10 x 10 grid for 1.0 mL cryo tubes	130x130x40 mm	Plastic	clear	case / 36 pcs.	3-636-45-0
Cryobox, 10 x 10 grid for 2.0 mL cryo tubes	130x130x50 mm	Plastic	clear	case / 36 pcs.	3-637-45-0
Cryobox, 10 x 10 grid for 4.0 mL cryo tubes	130x130x80 mm	Plastic	clear	case / 36 pcs.	3-638-45-0
Cryobox, 10 x 10 grid for 5.0 mL cryo tubes	130x130x90 mm	Plastic	clear	case / 36 pcs.	3-639-45-0

Cardboard Cryoboxes

Construction – Boxes and dividers are constructed from a temperature-resistant recyclable cardboard

Dimensions – 133x133x50 mm


Low cost – removable grid for low cost storage of your cryogenic samples

Cell dividers – Supplied as 81 (9x9) or 100 (10x10) grid options

Flexible storage options – removable inserts can be swapped for different purposes, or removed completely for larger vessels

Water-repellent finish – for durability

Recyclable cardboard

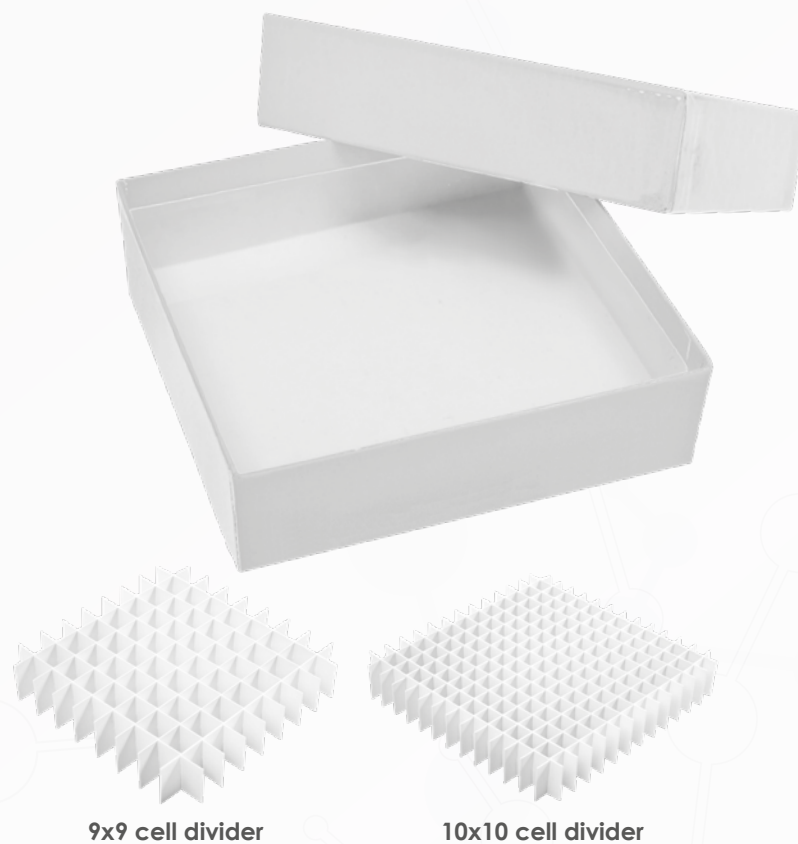

Ordering information

Cardboard Cryoboxes

Description	Dimensions	Material	Colour	Sales Unit	Cat. No.
Cryobox B9, includes 9x9 cell dividers	133x133x50 mm	Cardboard	white	case / 36 pcs.	3-640-36-9
Cryobox B10, includes 10x10 cell dividers	133x133x50 mm	Cardboard	white	case / 36 pcs.	3-641-36-9

For more colour and cell divider options please contact info@cappahn.com

Additional colour options 



ahn myTube | myPlate



Purification Systems

24

ahn myTube SC Spin Column System

Spin column system - consisting of 0.8 mL filter tube (round bottom design) equipped with corresponding filter (see filter options below) and 2.0 mL receiver tube

Elimination of contamination and loss of sample - high quality raw material meets purity standards, highly resistant to extreme mechanical strain and has an excellent chemical resistance

Centrifuge specifications: up to 25.000 x g
- fixed angle rotor up to 30.000 x g - swing out rotor



ahn myPlate FP Microfilterplate System

Microfilterplate - 8x12 matrix, 96 individual drip directors, 1.0 mL capacity, in accordance with SBS-Standard, equipped with corresponding filter (see filter options below) / receiver plate and vacuum manifold are separately available

Elimination of risk of sample contamination - manufactured from absolutely clear high-quality virgin polystyrene

Filter Options	SC	FP
Glass fibre (GF) filter characteristics - recommended for the purification of DNA and RNA / prefiltering of dirty solutions / rapid flow rate and high particle loading capacity	■	■
Polyethylene (PE) filter characteristics - suitable for particle-removing coarse filtration / can be used both as a prefilter or as support filter in combination with other filters / hydrophobic / very good chemical resistance	■	■

Ordering information



Description - spin column system

Description	Colour	Pack Type	Sales Unit	Cat. No.	
AHN myTube® SC 0.8 mL filter tube, GF-F2 filter, 2.0 mL receiver tube	①	clear	bag	2x100	3-320-20-0
AHN myTube® SC 0.8 mL filter tube, GF-N2 filter, 2.0 mL receiver tube	①	clear	bag	2x100	3-321-20-0
AHN myTube® SC 0.8 mL filter tube, PE 20-60 µm filter, 2.0 mL receiver tube	①	clear	bag	2x100	3-340-10-0

Description - individual components spin column system

Description	Colour	Pack Type	Sales Unit	Cat. No.	
AHN myTube® SC 0.8 mL filter tube w/o filter, w/filter ring	②	clear	bag	2x1000	3-331-50-0
AHN myTube® RT 1.5 mL receiver tube	③	clear	bag	10x1000	3-200-80-0
AHN myTube® RT 2.0 mL receiver tube	④	clear	bag	10x1000	3-202-80-0

Description - microplate system

Description	Colour	Pack Type	Sales Unit	Cat. No.	
AHN myPlate® FP 96-wells, 1.0 mL, GF-FF filter	⑤	clear	bag	5x5	3-210-25-0
AHN myPlate® FP 96-wells, 1.0 mL, GF/N 1.6 µm filter	⑤	clear	bag	5x5	3-210-25-2
AHN myPlate® FP 96-wells, 1.0 mL, PE 20-60 µm filter	⑤	clear	bag	5x5	3-225-25-0

Note: Further filter options are available on request. These include Nylon (recommended for filtering protein-free culture media), Cellulose Acetate (recommended for applications requiring low protein binding), Regenerated Cellulose (suitable for particle removal from solvents) and Polyvinylidene Fluoride (used in general filtration and sample preparation applications). Contact info@cappahn.com for ordering information.



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



ahn pipet4U pro

Manual Single and Multi Channel Pipettes

26

Comfortable and effortless pipetting - new user-friendly handle design with soft grip for comfortable operation to prevent hand stress / light weight / low plunger force minimises RSI (repetitive strain injury)

Easy identification - use the easily attachable colour-coded drops to identify your pipette according to your requirements

Durability - new tip cone demonstrates high mechanical strength and high resistance against chemicals

Calibration stability - even after multiple autoclaving no recalibration is required


Accurate and reliable pipetting - reverse pipetting option for accurate and precise aspirating and dispensing of viscous liquids / innovative soft grip prevents transfer of heat to internal components, ensuring accurate pipette operation / lock counter mechanism of variable volume pipettes prevents volume changes while pipetting / blow out system allows last drop dispensing

Sterile handling - fully autoclavable pipette (by steam at 121°C / 15 min)

Easy and quick maintenance and adjustment - with the included disassembly and calibration tool

Meets following standards: DIN EN ISO 8655-1
| DIN EN ISO 8655-2 | DIN EN ISO 8655-6
| EU Directive 98/79/EG

Compatibility list available on request

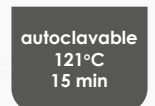

Colour Coding



Ergonomic Shape



Delivery Package

1 pc. Manual pipette	1 bag Sample Tips
1 pc. Calibration / Disassembling Tool	1 bag O-rings
1 pc. Drop Tool	1 pc. Instruction Manual
1 bag Colour Coding Drop	1 pc. Warranty Card
1 pc. Silicon Grease	1 pc. Calibration Report

ahn



Uthleber Weg 14 | 99734 Nordhausen | Germany

Ordering information

Manual single channel pipettes, fixed volume

Volume	Test volume	Increments	Inaccuracy ¹	Imprecision ²	Cat. No.
5 µL	5 µL	-	1.30%	1.20%	8-001-20-9
10 µL	10 µL	-	0.80%	0.80%	8-002-20-9
20 µL	20 µL	-	0.60%	0.50%	8-003-20-9
25 µL	25 µL	-	0.50%	0.30%	8-004-20-9
50 µL	50 µL	-	0.50%	0.30%	8-007-20-9
100 µL	100 µL	-	0.50%	0.30%	8-008-20-9
200 µL	200 µL	-	0.40%	0.20%	8-009-20-9
250 µL	250 µL	-	0.40%	0.20%	8-010-20-9
500 µL	500 µL	-	0.30%	0.20%	8-011-20-9
1000 µL	1000 µL	-	0.30%	0.20%	8-012-20-9
2 mL	2 mL	-	0.30%	0.15%	8-013-20-9
2.5 mL	2.5 mL	-	0.30%	0.15%	8-014-20-9
5 mL	5 mL	-	0.30%	0.15%	8-015-20-9
10 mL ³	10 mL	-	0.60%	0.20%	8-016-20-9

Manual single channel pipettes, variable volume

Volume	Test volume	Increments	Inaccuracy ¹	Imprecision ²	Cat. No.
0.1-2.5 µL	0.25 µL		12.00 %	6.00 %	8-100-20-9
	1.25 µL	0.01 µL	3.00 %	3.00 %	
	2.5 µL		2.50 %	1.60 %	
0.5-10 µL	1 µL		2.50%	1.50%	8-101-20-9
	5 µL	0.10 µL	2.00%	1.00%	
	10 µL		1.00%	0.80%	
2-20µL	2 µL		3.00%	2.00%	8-102-20-9
	10 µL	0.10 µL	1.50%	1.00%	
	20 µL		0.90%	0.40%	
5-50 µL	5 µL		2.00%	2.00%	8-103-20-9
	25 µL	0.50 µL	0.80%	0.40%	
	50 µL		0.60%	0.30%	
10-100 µL	10 µL		3.00%	1.50%	8-104-20-9
	50 µL	0.50 µL	1.00%	0.50%	
	100 µL		0.80%	0.15%	
20-200 µL	20 µL		2.00%	0.80%	8-105-20-9
	100 µL	1.00 µL	0.70%	0.30%	
	200 µL		0.60%	0.15%	
50-200 µL	50 µL		1.00%	0.40%	8-105-21-9
	100 µL	1.00 µL	0.80%	0.20%	
	200 µL		0.60%	0.15%	
100-1000 µL	100 µL		2.00%	0.70%	8-106-20-9
	500 µL	5.00 µL	1.00%	0.40%	
	1000 µL		0.60%	0.20%	
200-1000 µL	200 µL		0.90%	0.30%	8-106-21-9
	500 µL	5.00 µL	0.75%	0.25%	
	1000 µL		0.60%	0.20%	
0.5-5 mL	0.5 mL		2.00%	0.60%	8-107-20-9
	2.5 mL	0.05 mL	0.60%	0.30%	
	5 mL		0.50%	0.15%	
2-10 mL ³	2 mL		3.00%	0.60%	8-108-20-9
	5 mL	0.10 mL	1.20%	0.30%	
	10 mL		0.60%	0.20%	

¹ Systematic error

² Random error

³ Tip ejector function not available

AHN pipet4u pro Pipette Starter Kits

Description	Content	Cat. No.
Pipette Starter Kit	4 pipettes (0.5-10 µL/10-100 µL/100-1000 µL/2-10 mL), a 6-place Pipette Stand, 3 racks of Pipette Tips and 1 bag of Pipette Tips	8-SK1-99-9
Pipette Starter Kit	4 pipettes (0.1-2.5 µL/2-20 µL/20-200 µL/ 100-1000 µL), a 6-place Pipette Stand, 4 racks of Pipette Tips	8-SK2-99-9










Ordering information

Manual multi channel pipettes, variable volume

Volume	Channels	Test volume	Increments	Inaccuracy ¹	Imprecision ²	Cat. No.
0.5-10 µL	8	0.5 µL	0.10 µL	4.00%	4.00%	8-201-20-9
		5 µL		2.50%	2.50%	
		10 µL		1.50%	1.50%	
5-50 µL	8	5 µL	0.50 µL	3.00%	2.00%	8-202-20-9
		25 µL		1.50%	1.00%	
		50 µL		1.00%	0.70%	
10-100 µL	8	10 µL	1.00 µL	1.50%	0.75%	8-204-20-0
		50 µL		1.00%	0.50%	
		100 µL		1.00%	0.25%	
30-300 µL	8	30 µL	1.00 µL	1.50%	0.75%	8-203-20-9
		150 µL		1.00%	0.50%	
		300 µL		0.80%	0.25%	
0.5-10 µL	12	1 µL	0.10 µL	4.00%	4.00%	8-401-20-9
		5 µL		2.50%	2.50%	
		10 µL		1.50%	1.50%	
5-50 µL	12	5 µL	0.50 µL	3.00%	2.00%	8-402-20-9
		25 µL		1.50%	1.00%	
		50 µL		1.00%	0.70%	
10-100 µL	12	10 µL	1.00 µL	1.50%	0.75%	8-404-20-0
		50 µL		1.00%	0.50%	
		100 µL		1.00%	0.25%	
30-300 µL	12	30 µL	1.00 µL	1.50%	0.75%	8-403-20-9
		150 µL		1.00%	0.50%	
		300 µL		0.80%	0.25%	

¹ Systematic error ² Random error

Spare parts for pipet4u pro manual single channel pipettes

Set	Volume	Cat. No.
	0.1-2.5 µL	8-190-00-9
	0.5-10 µL	8-190-01-9
	2-20 µL	8-190-02-9
	5 µL	8-190-03-9
	10 µL	
	20 µL	
	25 µL	
	50 µL	
	5-50 µL	
	100 µL	8-190-04-9
	10-100 µL	
	200 µL	8-190-05-9
	20-200 µL	
	50-200 µL	
	250 µL	8-190-06-9
	500 µL	
	1000 µL	
	100-1000 µL	
	200-1000 µL	
	2 mL	8-190-07-9
	2.5 mL	
	5 mL	
	0.5-5 mL	
	10 mL	
	2-10 mL	8-190-08-9


ahn Pipette Stands

Single and Multi Channel Pipettes

Safe and effortless handling - comfortable and reliable pipette storage / pipette stands can be easily assembled and disassembled for cleaning purposes

Versions - available as 1 place / 4 places and 6 places pipette stands for single channel and multi channel pipettes and as fixed place and variable place pipette stands


Acrylic / Aluminium pipette stands

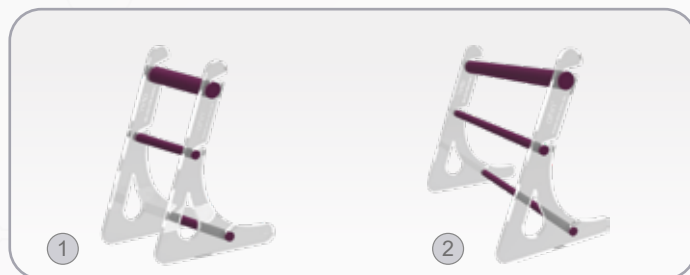

1

2

Acrylic pipette stands

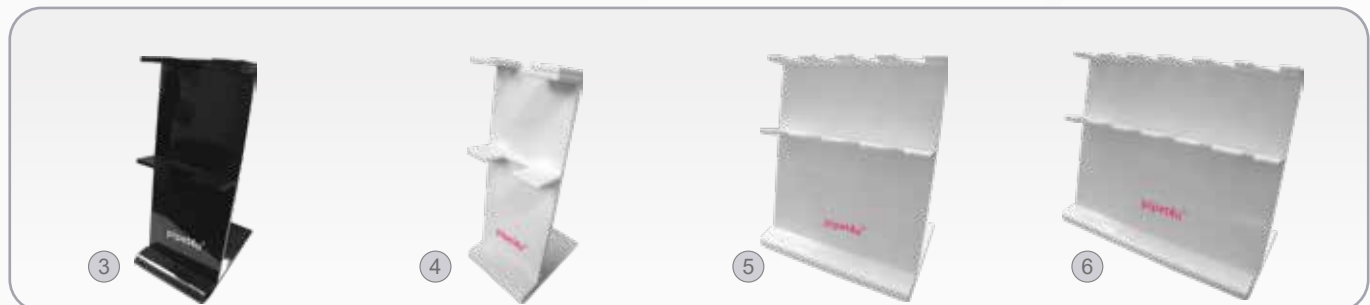

3

4

5

6



Ordering information

Description - Acrylic/Aluminium pipette stands		Colour	Sales Unit	Cat. No.
1 place pipette stand for manual single and multi channel pipettes variable place	①	clear/magenta	Carton/1 pc.	8-132-01-8
6 places pipette stand for manual single channel pipettes variable place	②	clear/magenta	Carton/1 pc.	8-132-00-8

Description - Acrylic pipette stands		Colour	Sales Unit	Cat. No.
1 place pipette stand for manual multi channel pipettes fixed place	③	black	Carton/1 pc.	8-130-02-7
1 place pipette stand for manual multi channel pipettes fixed place	④	white	Carton/1 pc.	8-130-02-9
4 places pipette stand for manual single channel pipettes fixed place	⑤	white	Carton/1 pc.	8-130-00-9
6 places pipette stand for manual single channel pipettes fixed place	⑥	white	Carton/1 pc.	8-131-00-9



ahn pipet4U Revolution

Electronic Single and Multi Channel Pipettes

30

Comfortable and effortless pipetting - ultra bright OLED display for extra clear visibility / ergonomic design for fatigue-free operation / comfortable pipetting dial for easy volume control and low strain pipetting (minimises RSI - repetitive strain injury)

Easy handling with one dial - control of all operating programs with the patented one dial system

6 different operating programs - pipetting / dispensing / reverse pipetting / mixing / autocalibration / speed control

Accurate and reliable pipetting - SCS-System (speed control system) controls piston movement giving highest precision and accuracy

Easy self-calibration - ESM-System (easy self-calibration management system) allows user calibration of the pipette, the user can calibrate their pipette for solutions of different viscosity, the system also controls permanently all functions and indicates errors

Sterile handling - autoclavable pipette cone (120°C / 15 min)

Meets following standards: DIN EN ISO 8655-1 | DIN EN ISO 8655-2 | DIN EN ISO 8655-6 | EU Directive 98/79/EG

Compatibility list available on request

Ordering information

Electronic single channel pipettes

Volume	Test volume	Increments	Inaccuracy ¹	Imprecision ²	Cat. No.
0.2-10 µL	0.5 µL	0.01 µL	5.00%	3.00%	8-500-00-0
	5 µL		2.50%	1.50%	
	10 µL		0.90%	0.50%	
5-100 µL	5 µL	0.10 µL	2.50%	1.50%	8-501-00-0
	50 µL		0.90%	0.30%	
	100 µL		0.40%	0.50%	
10-200 µL	10 µL	0.10 µL	0.90%	0.50%	8-502-00-0
	100 µL		0.40%	0.50%	
	200 µL		0.30%	0.10%	
50-1000 µL	100 µL	1.00 µL	0.40%	0.50%	8-503-00-0
	500 µL		0.40%	0.20%	
	1000 µL		0.30%	0.10%	
100-5000 µL	500 µL	5.00 µL	0.40%	0.50%	8-505-00-0
	2500 µL		0.10%	0.10%	
	5000 µL		0.10%	0.10%	

¹ Systematic error ² Random error

Delivery Package

1 pc. Electronic pipette	1 pc. Li-Ion-battery
1 pc. Sample bag with tips	1 pc. Charger (EU)
1 pc. Grease	1 pc. Instruction Manual
1 bag Screwdriver	1 pc. Calibration Report

UV
RESISTANT

CE IVD

2
YEARS
WARRANTY

autoclavable
120°C
15 min



Ordering information

Electronic multi channel pipettes

Volume	Channels	Test volume	Increments	Inaccuracy ¹	Imprecision ²	Cat. No.
0.2-10 µL	8	0.5 µL	0.01 µL	5.00%	3.00%	8-520-00-0
		5 µL		2.50%	1.50%	
		10 µL		0.90%	0.50%	
2.5-50 µL	8	5 µL	0.10 µL	2.50%	1.50%	8-521-00-0
		25 µL		0.90%	0.30%	
		50 µL		0.40%	0.50%	
10-100 µL	8	10 µL	0.10 µL	3.00 %	2.00 %	8-523-00-0
		50 µL		1.00 %	0.80 %	
		100 µL		0.80 %	0.30 %	
15-300 µL	8	20 µL	0.50 µL	0.90%	0.50%	8-522-00-0
		150 µL		0.40%	0.30%	
		300 µL		0.30%	0.10%	
0.2-10 µL	12	0.5 µL	0.01 µL	5.00%	3.00%	8-530-00-0
		5 µL		2.50%	1.50%	
		10 µL		0.90%	0.50%	
2.5-50 µL	12	5 µL	0.10 µL	2.50%	1.50%	8-531-00-0
		25 µL		0.90%	0.30%	
		50 µL		0.40%	0.50%	
5-100 µL	12	10 µL	0.10 µL	3.00 %	2.00 %	8-533-00-0
		50 µL		1.00 %	0.80 %	
		100 µL		0.80 %	0.30 %	
15-300 µL	12	20 µL	0.50 µL	0.90%	0.50%	8-532-00-0
		150 µL		0.40%	0.30%	
		300 µL		0.30%	0.10%	

Delivery Package

1 pc. Electronic pipette	1 pc. Li-Ion-battery
1 pc. Sample bag with tips	1 pc. Charger (EU)
1 pc. Grease	1 pc. Instruction Manual
1 bag Screwdriver	1 pc. Calibration Report

¹ Systematic error ² Random error



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.dhn-bio.de



ahn pipet4u Oasis

Pipette Controller

32

Comfortable and effortless pipetting - light weight pipette design / smooth push buttons for precise pipetting control / ergonomic design for easy operation / cordless working (holds charge up to 8 hours) / quiet high-powered motor

Contamination control - replaceable membrane filter to minimise risk of contamination

Speed setting - adjustment of pipetting speed on sample viscosity with 3 selectable speeds

Sterile handling - autoclavable pipette cone (120°C / 15 min)



UV RESISTANT	CE IVD	2 YEARS WARRANTY	autoclavable 120°C 15 min
---------------------	---------------	-------------------------	----------------------------------

Ordering information


Description	Sales Unit	Cat. No.
pipet4u® Oasis Pipette Controller 0.1-100 mL	Carton/1 pc.	8-300-00-9

Description - Spare Parts	Sales Unit	Cat. No.
0.20 µm PTFE membrane filter	Bag/5 pcs.	8-300-03-0
0.45 µm PTFE membrane filter	Bag/5 pcs.	8-300-04-0
Silicone adapter	Bag/1 pc.	8-300-01-0
Pipette cone	Bag/1 pc.	8-300-02-8

Delivery Package
1 pc. Pipette Controller
2 pcs. PTFE Membrane filters (0.2 µm / 0.45 µm)
1 pc. Charger (EU)
2 pcs. NiMH batteries
1 pc. Wall mount bracket
1 pc. Desktop stand
1 pc. Instruction manual
1 pc. Warranty card



Autoclavable tip cone

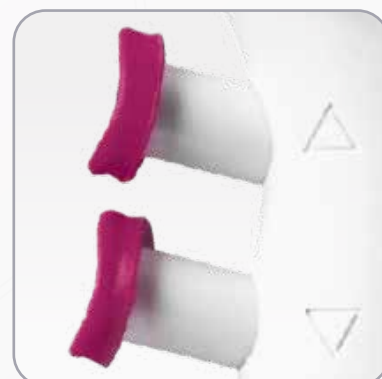

Ergonomic push buttons

ahn myDpette

Bottle Top Dispenser

Unique piston mechanism - cleaning of piston and cylinder without affecting the calibration

Dual graduation marks - readability of the volume from both sides

Springless valve design - leads to smooth functioning and high chemical resistance

Eco-friendly - while discharging the tubes in the re-circulation mode, in fast priming, the reagent is being redirected back into the bottle, eliminating any waste of precious reagents, ensuring safety and savings

Easy handling - discharge tube with 360° rotation, makes the bottle label visible all the time / telescopic filling tube compatible with different sizes of reagent bottles / ergonomic grip

Easy to calibrate - complies with ISO 8655-5 standard



33

Ordering information

Description	Sales Unit	Cat. No.
AHN myDpette Bottle Top Dispenser 0.25-2.5 mL	Carton/1 pc.	8-700-00-0
AHN myDpette Bottle Top Dispenser 0.5-5.0 mL	Carton/1 pc.	8-701-00-0
AHN myDpette Bottle Top Dispenser 1.0-10.0 mL	Carton/1 pc.	8-702-00-0
AHN myDpette Bottle Top Dispenser 2.5-25.0 mL	Carton/1 pc.	8-703-00-0
AHN myDpette Bottle Top Dispenser 5.0-50.0 mL	Carton/1 pc.	8-704-00-0

2
YEARS
WARRANTY

autoclavable
120°C
20 min

Delivery Package

1 pc. Bottle top dispenser
4 pcs. Adapters (28, 38, 40 and 45 mm screw cap)
1 pc. Telescopic filling tube
1 pc. Re-circulation tube
1 pc. Calibration tool
1 pc. User manual
1 pc. Certificate of conformity
1 pc. Calibration report
1 pc. Warranty card



Double scale



Recirculation valve



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



ahn myLab MC-01

Microcentrifuge 6000 rpm

34

Usage - separation of particles from suspension or macromolecules from solutions

Capacity - 8-place closed rotor for 1.5 mL / 2.0 mL micro tubes

Maximum Speed - 6000 rpm / 2000 x g

Multiple tube options - adaptors for 0.2 mL / 0.5 mL tubes available

Unique closed rotor design - ensures better airflow through low-friction, stability for temperature sensitive samples through low heat generation and ultra-silent, vibration free operation


Secure operation - electronic brake for immediate stop on lid opening for increased user safety

Durability - robust construction / solid hinge and switch design

Quick and comfortable operation - possibility of quick runs by simply closing and opening the lid / digital display

Small footprint saves valuable bench space - suitable for use in hoods and cold rooms

Speed accuracy +/- 100 rpm

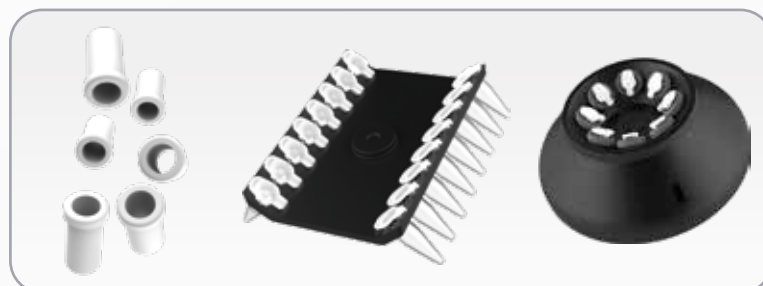

Ordering information

Description	Cat. No.
AHN myLab® MC-01 Microcentrifuge 6000 rpm	7-000-01-0

Specifications	
Noise level	< 55 db
Quick spin feature	yes, lid activated on/off
Power consumption	10 W
Dimensions (WxDxH)	162x157x155 mm

Delivery Package
1 pc. Microcentrifuge 6000 rpm
1 pc. 1.5 / 2.0 mL closed rotor
1 pc. PCR strip rotor
1 pc. T-Allen wrench
8 pcs. Individual adaptor for 0.2 mL tubes
8 pcs. Individual adaptor for 0.5 mL tubes
1 pc. External power adaptor
1 pc. Instruction manual
1 pc. Warranty card
1 pc. Declaration of conformity

ahn myLab MC-02

High Speed Microcentrifuge 15000 rpm

Usage - designed for blood centrifugation and other clinical applications

Capacity - 12-place closed rotor for 1.5 mL / 2.0 mL micro tubes

Maximum Speed - up to 15000 rpm / 15596 x g

Multiple tube options - adaptors for 0.2 mL / 0.5 mL tubes available

Unique closed rotor design - ensures better airflow through low-friction and stability for temperature sensitive samples through low heat generation

Increased user safety and comfort - automatic safety brake system preventing the lid from opening upon operation / imbalance detection / rotor with screw on lid to prevent contamination

Small footprint - optimizing bench top space

Multiple applications - suitable for use in hoods or cold rooms

Robust design - brushless DC motor enables maintenance-free and long-lasting operation, perfect for extended runs

Comfortable use - user-friendly interface / extremely low noise level / large digital display



Ordering information

Description	Cat. No.
AHN myLab® MC-02 Microcentrifuge 15000 rpm	7-001-01-0

Optional attachments	Cat. No.
PCR strip rotor	CR-15-PCR

Specifications	
Noise level	< 60 db
Imbalance detection	yes, auto cut off
Timer setting	1 to 99 mins / infinite mode
Dimensions (WxDxH)	260x244x205 mm



PCR strip rotor

Delivery Package	
1 pc.	Microcentrifuge 15000 rpm
1 pc.	Power supply adaptor
1 pc.	12 slot rotor (pre-installed)
1 pc.	T-Allen key
12 pc.	Reduction adaptors for 0.2 / 0.5 mL tube
1 pc.	Instruction manual
1 pc.	Warranty card



ahn myLab CLC-01

Clinical Centrifuge 4000 rpm

36

Usage - designed for blood centrifugation and other clinical applications

Multiple rotor and tube options - 6x10 mL (swing out rotor) / 8x15 mL (fixed angle rotor) / 16x10 mL (fixed angle rotor) / with each rotor, the appropriate tube adaptors, sleeves, spacers are included - 10 mL, 9 mL, 8 mL, 6 mL, 4 mL, 3 mL and 2 mL

Maximum Speed - up to 4000 rpm / 2270 x g varies per rotor

Precise operation - speed interval of 10 rpm starting from 500 rpm

Memory feature - can save up to 99 user-defined programs (protocols)




Increased user safety and comfort - automatic safety brake system preventing the lid from opening upon operation / imbalance detection

Robust design - brushless DC motor enables maintenance-free and long-lasting operation, perfect for extended runs

Comfortable use - last run memory function / user-friendly interface for speed and time setting / quick spin feature



Ordering information

Specifications					
Rotor Type: Swing out		Rotor Type: Fixed angle		Rotor Type: Fixed angle	
Capacity: 6x10 mL		Capacity: 8x15 mL		Capacity: 16x10 mL	
Max RCF: 2270 x g		Max RCF: 1950 x g		Max RCF: 2180 x g	
Cat. No.: CRC-4-6S		Cat. No.: CRC-4-8		Cat. No.: CRC-4-16	

Description	Cat. No.
AHN myLab® CLC-01 Clinical Centrifuge 4000 rpm	7-010-00-0

Specifications	
Imbalance detection	yes, auto cut off
Timer setting	1-999 mins / infinite mode
Weight	7.7 kg (without rotor)
Dimensions (WxDxH)	355x415x173 mm

Delivery Package
1 pc. Clinical centrifuge 4000 rpm
1 pc. Power supply adaptor
1 pc. T-Allen Key
1 pc. Swing out rotor
2 x 6 pcs. Tube holders
1 pc. Instruction manual
1 pc. Warranty card

AHN myLab CLC-01 Clinical Centrifuge 4000 rpm, is designed for processing blood and other clinical samples in applications like PRP, medical practitioners, pathological laboratories, as well as for routine hospital centrifuging test. In particular it is used for processing and analyzing samples from human body in in vitro diagnostic applications to ensure that in vitro diagnostic device can be used according to Directive 98/79/EC. Therefore AHN myLab CLC-01 Clinical Centrifuge 4000 rpm, and its accessories being in-vitro diagnostic accessories are themselves in-vitro diagnostic medical devices within the meaning of directive 98/79/EC of the European Parliament and of the Council of 27 October 1998 on in-vitro diagnostic medical devices.

ahn myLab HC-01

Hematocrit Centrifuge 12000 rpm

Usage - determination of volume fractions of erythrocytes (red blood cells) in blood and also for separation of blood and solutions

Capacity - 24 capillary tubes of 75 / 45 mm length

Maximum Speed - up to 12000 rpm / 14119 x g

Precise operation - speed interval of 100 rpm starting from 500 rpm

Increased user safety and comfort - automatic safety brake system preventing the lid from opening upon operation / imbalance detection / HTC rotor with screw on lid to prevent contamination

User comfort - big user-friendly display / last run memory function / preset feature where the user can set up their favorite mode and change the speed (rpm) and time

Robust design - brushless DC motor enables maintenance-free and long-lasting operation, perfect for extended runs

Small footprint - optimizing bench top space



37



24 slot capillary rotor

Ordering information

Description	Cat. No.
AHN myLab® HC-01 Hematocrit Centrifuge 12000 rpm	7-020-00-0

Specifications	
Imbalance detection	yes, auto cut off
Timer setting	up to 12 mins
Weight	5.0 kg (without Rotor)
Dimensions (WxDxH)	301x363x105 mm

Delivery Package
1 pc. Hematocrit Centrifuge 12000 rpm
1 pc. Power supply adaptor
1 pc. 24-slot capillary rotor (pre-installed)
1 pc. Blood reader
1 pc. T-Allen wrench
1 pc. Instruction manual
1 pc. Warranty card



+49(0)3631/65242-0



+49(0)3631/65242-90



info@cappahn.com



www.ahn-bio.de



ahn myLab UC-01

Universal Centrifuge 4500 rpm

38

Wide range of rotor options guarantee the most versatile scope of application

User-friendly interface - program the device according to your specific needs.

Brushless DC motor technology sustains durability and high performance over many years of use

Adjustable speeds - range from 500 to 4.500 rpm

Multiple rotor options - with both swing-out and fixed angle rotors

Large digital display - visualises multiple parameters in a clear, comfortable way

9 acceleration and 9 deceleration modes - enabling you to start and stop the device without disturbing fragile samples

Adjustable timer - can be set from 1 to 999 minutes, or the infinite centrifugation mode which has no time limitation

Excellent user safety - thanks to the electronic imbalance detection feature, which is activated by uneven loading of samples

Lid Lock Safety Break - ensures the lid remains closed until the rotor has fully stopped. The Emergency Lid Release only allows you to open the device once electricity supply has been removed

Short spin feature - ensures the ultimate performance for ultra-short centrifugation and pelleting. The stainless steel chamber rotor is fully corrosion-free and easy to clean

Intelligently designed airflow system prevents increase in sample temperature.

Small footprint - uses minimal bench space while offering powerful and efficient centrifugation and maintaining a high performance



Specifications

Speed Range	500 - 4.500 rpm / (g) varies per rotor type
Motor Type	Brushless DC motor (BLDC)
Remote Operation	Yes, via USB
Display	Yes, large LED Display
Timer Setting	from 1 to 999 mins or infinite mode
Programmable	Yes, up to 99 programs
Imbalance Detection	Yes
Rotors (optional)	Multiple rotor options (swing-out and fixed angle rotors)
Dimensions (WxDxH)	595x436x311 mm
Weight	23 kg (without rotor)

CR-43-100

4 x 100 mL swing rotor
with Bio Safety Cap,
RCF 3164 x g

**CR-43-50**

6 x 50 mL swing rotor
with
stainless steel tubes,
RCF 3530 x g

**CR-43-15**

16 x 15 mL swing rotor
with stainless steel
tubes, RCF 3485 x g

**CR-43-PL**

2 x microplates swing
rotor, RCF 2355 x g

CR-43-50F

6 x 50 mL fixed angle
rotor, RCF 2644 x g

**CR-43-100F**

4 x 100 mL fixed angle
rotor, RCF 2559 x g

**CR-43-15F-24**


24 x 15 mL fixed angle
rotor, RCF 2852 x g

**CR-43-15F-16**

16 x 15 mL fixed angle
rotor, RCF 3144 x g

Ordering information

Description	Cat. No.
AHN myLab® UC-01 Universal Centrifuge 4500 rpm	7-030-00-0
4x100mL swing rotor with Bio Safety Cap, RCF 3164 x g	CR-43-100
6x50mL swing rotor with stainless steel tubes	CR-43-50
16x15mL swing rotor with stainless steel tubes, RCF 3485 x g	CR-43-15
2 x Microplates Swing rotor, RCF 2355 x g	CR-43-PL
4x100mL Fix angle rotor, RCF 2559 x g	CR-43-100F
6x50mL Fix angle rotor, RCF 2644 x g	CR-43-50F
24x15mL Fix angle rotor, RCF 2852 x g	CR-43-15F-24
16x15mL Fix angle rotor, RCF 3144 x g	CR-43-15F-16
30x15mL Fix angle rotor, RCF 2729 x g	CR-43-15F-30
32x10 mL swing rotor, RCF 3164 x g (includes 32x10 mL stainless steel tubes + 32x10 mL rubber adapters for tube dia.13x length 100 mm + 32 x adapters for tube dia. 13x length 75 mm	CR-43-10

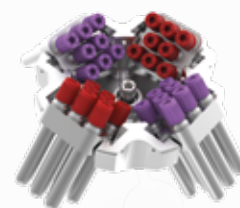**CR-43-10**

32 x 10 mL swing rotor,
RCF 3164 x g

**CR-43-15F-30**

30 x 15 mL fixed angle
rotor, RCF 2729 x g

Important note: Each rotor has a range of optional tube adaptors.
Please contact info@cappahn.com to get the full overview of all adaptors that are available.



+49(0)3631/65242-0



+49(0)3631/65242-90

info@cappahn.comwww.ahn-bio.de

ahn myLab VT-01

Mini Vortex Mixer 4500 rpm

40

Usage – quick mixing, re-suspension of samples in individual test tubes, up to a diameter of 30 mm

Maximum Speed – adjustable up to 4500 rpm

Robust design – brushless DC motor enables maintenance-free and long-lasting operation

Stable operation – through aluminium die cast base, suction cups and inbuilt counter balance system

Very small footprint saves valuable bench space - suitable for use in hoods and cold rooms



Soft silicone cup



Speed knob



Ordering information

Description	Cat. No.
AHN myLab® VT-01 Mini Vortex Mixer 4500 rpm	7-100-00-0

Specifications	
Max. tube diameter	30 mm
Orbital diameter	4.5 mm
Max. loading capacity	500 g
Press activation	yes
Protection class	IP 43
Dimensions (WxDxH)	100x100x67 mm
Weight	0.55 kg

Delivery Package
1 pc. Mini Vortex Mixer 4500 rpm
1 pc. Power supply adaptor
1 pc. Instruction manual
1 pc. Warranty card

ahn myLab VT-02 / 03

Vortex Mixer 2800 / 4800 rpm

Usage – comfortable everyday mixing operations

Maximum Speed – adjustable up to 2800 rpm (7-110-00-0) / 4200 rpm (7-120-00-0)

Comfortable use – selection of different modes (touch, stand-by and continuous), depending on the desired use / smooth and silent operation

Stable operation – inbuilt counter balance system

Digital display (VT-03) - adjustable speed setting with an interval of 10 rpm and time setting from 0 - 999 minutes / displays speed and time alternatively at every 5 seconds

Programmable pulse mode - ON time and OFF time can be individually modified

Multiple applications – optional attachments for micro plate, micro tubes and flask sticker for permissible speed for various attachments

Robust design – brushless DC motor enables maintenance-free and long-lasting operation

Small footprint saves valuable bench space – for use in hoods in the temperature range of 5 to 40°C



41



Optional attachments*

Foam for microtubes

Foam for micro plate

Rubber holder for flask

Ordering information

Description	Cat. No.
AHN myLab® VT-02 Vortex Mixer 2800 rpm	7-110-00-0
AHN myLab® VT-03 Digital Vortex Mixer 4200 rpm	7-120-00-0

Delivery Package
1 pc. Vortex Mixer
1 pc. Power supply adaptor
1 pc. Standard cup head (pre-installed)
1 pc. Instruction manual
1 pc. Warranty card

Specifications	
Orbital diameter	4 mm
Different modes	touch, continuous
Speed accuracy	± 5.0%
Max. loading capacity	500 g
Protection class	IP21
Power supply	220, 240 VAC, 50 Hz
Dimensions (WxDxH)	205x136x138.5 mm
Weight	3.0 kg

Optional attachments*	Cat. No.
Universal platform	7-110-01-0
Foam for micro tubes*	7-110-02-0
Foam for micro plates*	7-110-03-0
Rubber holder for flask*	7-110-04-0

* only applicable in combination with the universal platform 7-110-01-0



ahn myLab SM-01

Single Plate Magnetic Stirrer 1500 rpm

42

Usage - in laboratories, pharmacies and schools, mixing liquid industry


Maximum speed - from 15 to 1500 rpm with an interval of 50 rpm / software-driven gradual acceleration

Time setting - ranging from 1 - 99 minutes / continuous operation feature

Improved mixing - pulse mode where the mixing direction changes every 30 seconds

Maintenance-free design - water spill, corrosion and chemically resistant / includes no moving parts / high quality stirring bar

User-friendly interface and design - digital display showing the real-time speed and the remaining time / last run memory function / extremely thin (only 16.0 mm thick) / flat surface enables easy wiping and cleaning

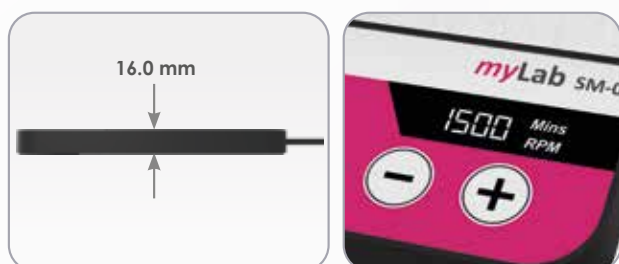

Ordering information

Description		Cat. No.	Specifications	
AHN myLab® SM-01 Single Plate Magnetic Stirrer 1500 rpm		7-300-00-0	Max. stirring quantity	800 mL
			Recommended stir bar length	25 mm
Optional accessories		Cat. No.	Protection class	IP65
Magnetic stirring bar - PTFE coated, 25 x 8 mm		CRS-MAG	Dimensions (WxDxH)	20x164x16 mm
			Weight	400 g
Delivery Package				
1 pc. Magnetic Stirrer 1500 rpm				
1 pc. Stir bar (PTFE coated, 25 x 8 mm)				
1 pc. Instruction manual				
1 pc. Warranty card				

ahn myLab SM-05/10/15

Heated Magnetic Stirrer 1500 rpm

Usage - chemical and pharma industry for mixing and / or heating liquids

Maximum Speed - adjustable speed ranging from 300 to 1500 rpm with an interval of 10 rpm

Multiple capacity options - 5, 10 or 15 flasks

Optimal heating process - heating up to 120°C with homogenous temperature distribution across the plate

Improved mixing - programmable pulse mode where the user can set the time from 30 to 99 seconds for changing the mixing direction

Time setting - ranging from 1 - 999 minutes / continuous operation feature

Operation-friendly design

- rubber feet for optimal stability on the bench / waterproof surface which is easy to clean and protects the equipment from accidental spills



Ordering information

Description	Cat. No.
AHN myLab® SM-05 5-plate Heated Magnetic Stirrer 1500 rpm	7-305-00-0
AHN myLab® SM-10 10-plate Heated Magnetic Stirrer 1500 rpm	7-310-00-0
AHN myLab® SM-15 15-plate Heated Magnetic Stirrer 1500 rpm	7-315-00-0

Optional accessories	Cat. No.
Magnetic stirring bar - PTFE coated, 25 x 8 mm	CRS-MAG

Specifications			
Number of stirring positions	5	10	15
Max. stirring quantity	400 mL per position		
Recommended stir bar length	25 mm		
Protection class	IP42		
Dimensions (WxDxH)	136x562x75 mm	196x565x75 mm	282x568x76 mm
Weight	3.0 kg	5.0 kg	7.0 kg

Delivery Package
1 pc. Heated Magnetic Stirrer 1500 rpm
1 pc. Power cord
1 pc. Reference stirrer bar
1 pc. Instruction manual
1 pc. Warranty card



ahn myLab TR-01

Tube Roller 30 rpm



Usage - mixing of blood samples and viscous substances

Capacity - 9 x 5 mL tubes

Fixed speed - 30 rpm with an accuracy of ± 2 rpm

Optimal mixing - gentle rocking and rolling movements at a tilt angle of ± 3 degree

Easy maintenance and adjustment - tool free removal and installation of rollers for effortless cleaning and replacing / easy adjustment of spacing to use with different tube diameter

Robust and space-saving design - robust construction for stable operation / small footprint requiring minimum of bench space



Tool-free removal and installation of rollers

Ordering information

Description	Cat. No.
AHN myLab® TR-01 Tube Roller 30 rpm	7-200-00-0

Specifications	
Run time	Continuous / infinite
Noise level	Noiseless operation
Roller size	30x330 mm
Max. loading capacity	1.5 kg
Protection class	IP 21
Dimensions (WxDxH)	475x201x80 mm
Weight	2.0 kg

Delivery Package
1 pc. Tube Roller 30 rpm
1 pc. Power supply adaptor
1 pc. Instruction manual
1 pc. Warranty card

ahn myLab TR-02

Tube Rotator 30 rpm

Usage - ideal for processing of blood samples

Capacity - 16 tubes (4 ESR Tubes and 12 blood collection tubes) with the maximum operational weight capacity of 1.5 kgs

Fixed speed - 30 rpm with an accuracy of ± 2 rpm

Optimal mixing - thanks to gentle circular rotations on the fixed plate at the angle of 38 degrees

Robust and space-saving design - robust metal construction with rubber feet for stable operation / compact design requires minimum of bench space

Easy maintenance and adjustment - the disc plate can be easily removed for cleaning and replacing



45

Ordering information

Description	Cat. No.
AHN myLab® TR-02 Tube Rotator 30 rpm	7-210-00-0

Specifications	
Run time	Continuous / infinite
Noise level	Noiseless operation
Max. loading capacity	1.5 kg
Protection class	IP 21
Dimensions (WxDxH)	278x183x252 mm
Weight	1.7 kg

Delivery Package
1 pc. Tube Rotator 30 rpm
1 pc. Power supply adaptor
1 pc. Instruction manual
1 pc. Warranty card



ahn myLab TR-03

Tube Rotator 30 rpm

46

Usage - ideal for processing of blood samples, separating and mixing of liquid and powdered samples

Capacity - 12 x 15 mL and 4 x 50 mL tubes or 22 x 15 mL tubes with the maximum operational weight capacity of 1.0 kg

Fixed speed - 30 rpm with an accuracy of ± 2 rpm

Optimal mixing - circular motion allows gentle but effective mixing

Robust and space-saving design - strong metal construction with rubber feet for stable operation / compact design requires minimum of bench space

Easy maintenance and adjustment - plates can be easily removed and installed for cleaning and replacing




Ordering information

Description	Cat. No.
AHN myLab® TR-03 Tube Rotator 30 rpm	7-220-00-0
Plate 22 x 15 mL (included with TR-03)	7-220-02-0

Optional attachments	Cat. No.
Plate 12 x 15 mL + 4 x 50 mL	7-220-01-0

Specifications	
Run time	Continuous / infinite
Noise level	Noiseless operation
Max. loading capacity	1.0 kg
Protection class	IP 21
Dimensions (WxDxH)	346x185x145 mm
Weight	2.5 kg

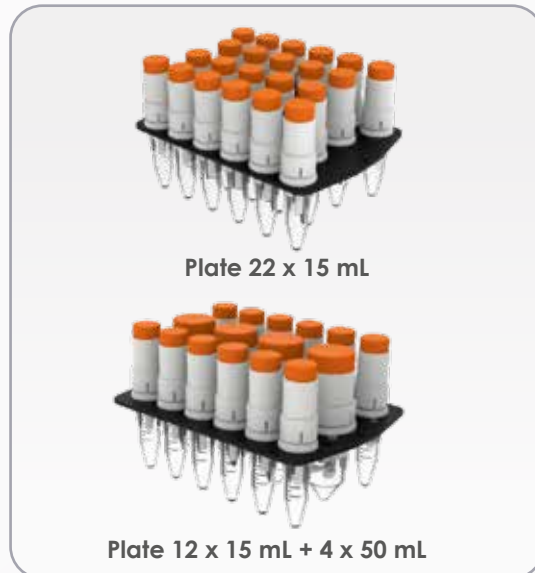

Delivery Package
1 pc. Tube Rotator 30 rpm
1 pc. Power supply adaptor
1 pc. Instruction manual
1 pc. Warranty card
1 pc. Plate 22 x 15 mL



+49(0)3631/65242-0



+49(0)3631/65242-90




info@cappahn.com



www.dhn-bio.de



Distributed by:

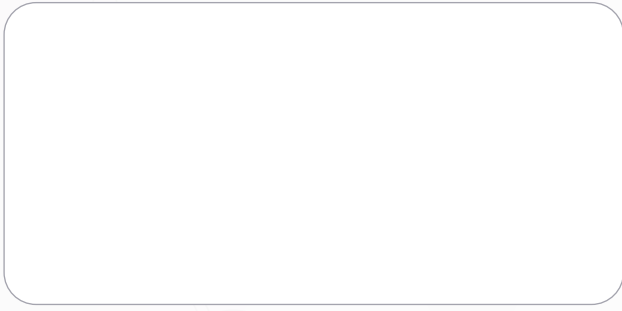

ahn

AHN Biotechnologie GmbH

Uthleber Weg 14
99734 Nordhausen
Germany

P: +49(0)3631/65242-0

F: +49(0)3631/65242-90

E: info@cappahn.com

www.ahn-bio.de

